


2025.08.19. <DNA Gel Excision>

- Agarose gel (LE agarose 2.25 g + 1X TAE Buffer 150 ml + SYBR green 4 ul)
- Electrophoresis at 100 V for 30 min.
- PCR product = 16 ul, 1kb plus DNA Ladder = 10 ul

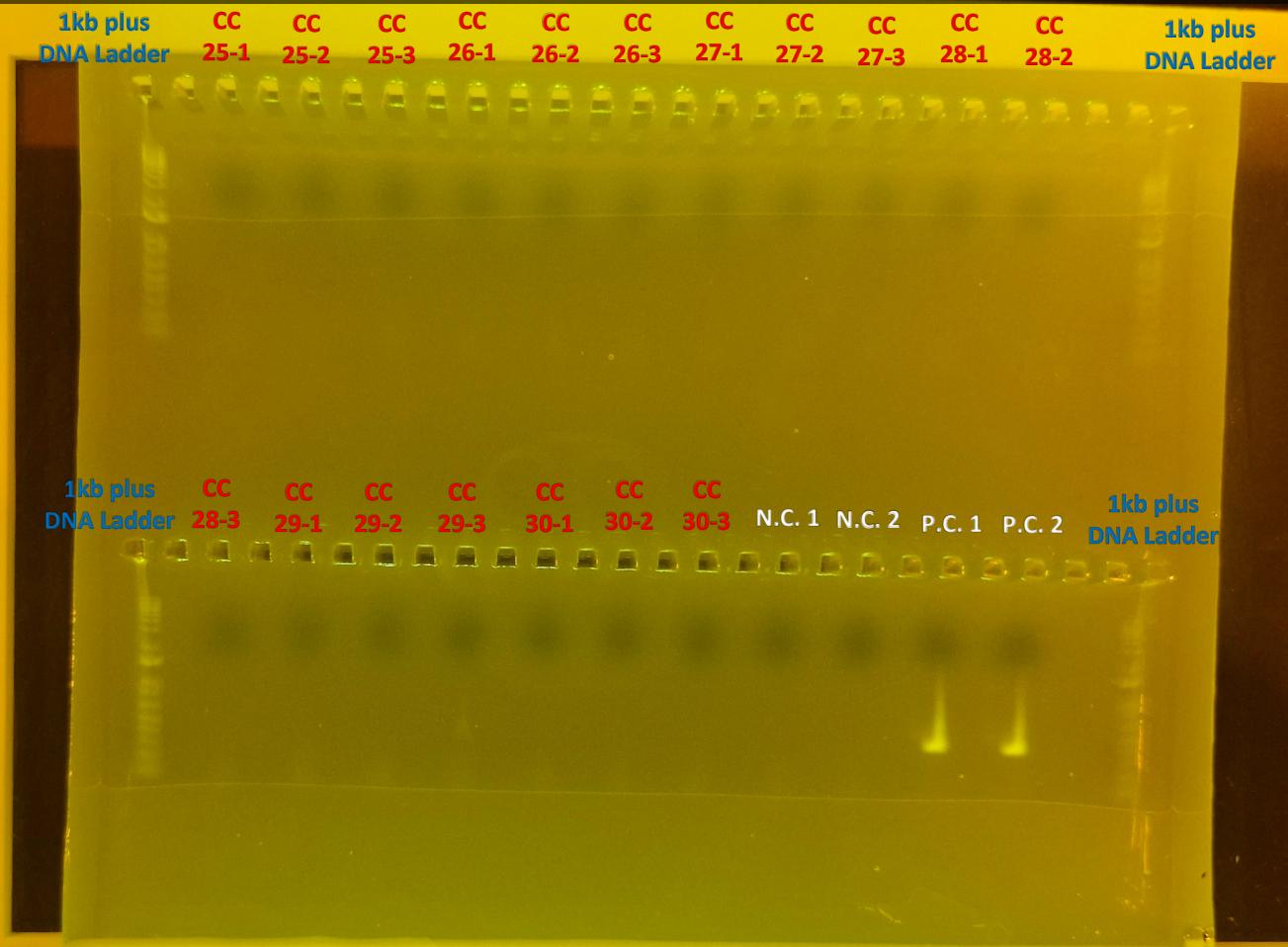

2025.08.20. <DNA Gel Excision>

- Agarose gel (LE agarose 2.25 g + 1X TAE Buffer 150 ml + SYBR green 4 ul)
- Electrophoresis at 100 V for 30 min.
- PCR product = 16 ul, 1kb plus DNA Ladder = 10 ul

1kb plus DNA Ladder CC 31-1 CC 31-2 CC 31-3 CC 32-1 CC 32-2 CC 32-3 CC 33-1 CC 33-2 CC 33-3 CC 34-1 CC 34-2 1kb plus DNA Ladder

1kb plus DNA Ladder CC 34-3 CC 35-1 CC 35-2 CC 35-3 CC 36-1 CC 36-2 CC 36-3 N.C. 1&2 P.C. 1&2 1kb plus DNA Ladder