


2025.07.15. <DNA Gel Excision>

- Agarose gel (LE agarose 2.25 g + 1X TAE Buffer 150 ml + SYBR green 4 ul)
- Electrophoresis at 100 V for 20 min.
- PCR product = 16 ul, 1kb plus DNA Ladder = 10 ul

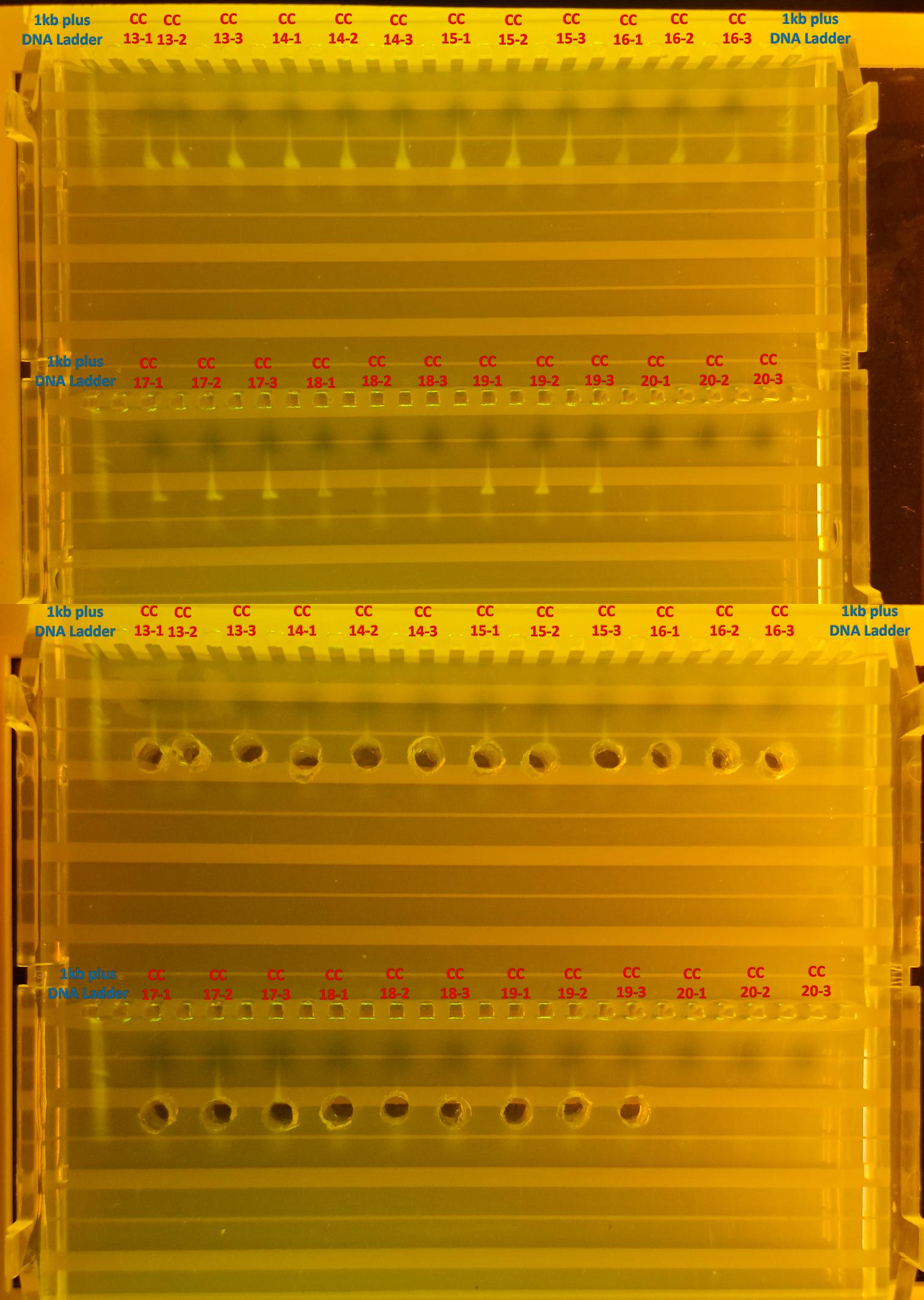

2025.07.15. <DNA Gel Excision>

- Agarose gel (LE agarose 2.25 g + 1X TAE Buffer 150 ml + SYBR green 4 ul)
- Electrophoresis at 100 V for 20 min.
- PCR product = 16 ul, 1kb plus DNA Ladder = 10 ul

